

COTTON PRODUCTS



To meet your need for perfection

... all things cotton

## roeko

- Reliable high manufacturing standard – qualityfocussed and fully automated
- Pharmaceutical quality raw materials – tested for purity, foreign fibres, absorbency and surfactant residues
- Versatile and well thought out – adapted to the needs of dentists

**Cotton products Made in Germany**

**roeko** – one name – long-standing experience – extensive technical know-how – utmost attention to quality

The quality of cotton products definitely pays dividends if one does the sums correctly.

The reliable high manufacturing standard of ROEKO cotton products contributes to lower consumption,

trouble-free working and hence to truly perfect results.



### LUNA GOLD

#### See and feel the difference

LUNA GOLD is encased in skin-friendly mesh and can be removed without prior wetting – this prevents mucosal reaction and saves time. The gold stripe clearly differentiates LUNA GOLD from conventional dental rolls.

#### For regular use we recommend


**LUNA** – long-lasting absorbency and exceptional dimensional stability

**COTTON-ROLL** – particularly soft with rapid absorption.

ROEKO dental rolls are bleached chlorine-free.

### ROLL-O-MAT

The system is both simple and ingenious: slight pressure on the lid is sufficient to dispense a single cotton roll and keep it positioned reliably on the dispenser tip. The ROLL-O-MAT is available as plastic dispenser or in a stainless steel version.

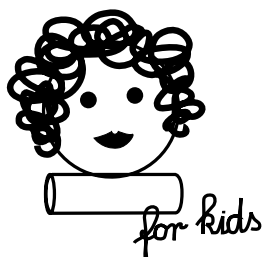


### LUNA FOR KIDS

#### Specially designed for the very small patients

LUNA FOR KIDS dental rolls are particularly short and yet dimensionally stable like their larger Luna counterparts. This leaves plenty of room to work, even in the smallest of mouths. ROEKO dental rolls are bleached chlorine-free.

### LUNAMAT

The LUNAMAT is equipped with a drawer which dispenses 2-3 cotton rolls when pulled out.

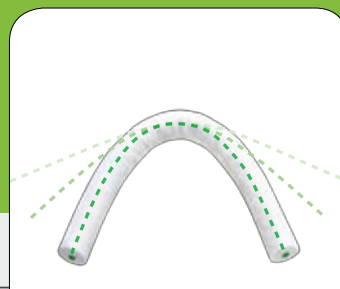

### PAROTISROLL

#### Multipurpose cotton roll with spring action

The spring rod inside the cotton roll provides a spring action and makes the PAROTISROLL self-fixating. When used buccally, lingually or labially, it keeps the mouth open, holds down the tongue and retracts the lips. When positioned appropriately, it also prevents the flow of saliva through pressure on the parotid gland.

### PAROTISROLL ANATOMIC

By Dr. Hannelore Oberschachtsiek  
PAROTISROLL anatomic is the Parotisroll with flattened ends. The flattened ends facilitate placement of the cotton roll and reduce the foreign body feeling in the mouth.

### PAROTIMAT

For convenient individual dispensing and hygienic storage of 50 PAROTISROLLS.



### DENTAL ROLL DISPENSERS

ROEKO dental roll dispensers keep cotton rolls within reach during treatment while at the same time protecting them hygienically from spray and dust.



## COTTON PELLETS

### Uniform in shape and size


The fully automated production of the ROEKO cotton pellets guarantees their uniform shape and size and thus, reliability. They are the ideal aid for drying cavities or applying medications.

#000	#00	#0	#1	#2
Ø 3 mm	Ø 4 mm	Ø 4,8 mm	Ø 6,3 mm	Ø 7,6 mm

## SOLOMAT

### The elegant dispenser

In this stainless steel dispenser, the cotton pellets are spring-loaded to the dispensing surface. The scratch- and abrasion-resistant surface prevents bacteria from adhering and is easy to clean and disinfect.

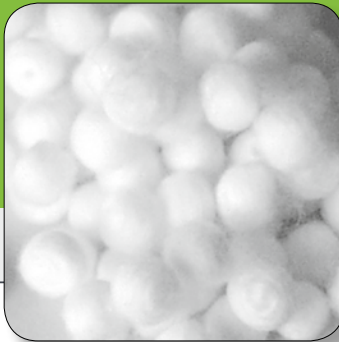

## SOLOMAT-N

### For individual dispensing of cotton pellets or foam pellets

Autoclavable plastic dispenser with refill system. The cotton pellets are moved manually to the dispensing surface.



## COMPRECAP

COMPRECAP compression caps are a simple and intelligent aid for preparing impression taking. They stop bleeding – by compression and dry tooth stumps.

## COMPRECAP ANATOMIC

COMPRECAP ANATOMIC are ideally adapted to the anatomy of the dental arch through semi-circular cut-outs.

## SOLOMAT-CLICK

### The smart dispensing system

SOLOMAT-CLICK consists of a dispensing grille, cover and refill. These 3 parts are clicked together when repurchasing a refill. The stainless steel base gives the dispenser the necessary stability.



## COTTON PELLETS DISPENSER


The ROEKO cotton pellet dispensers keep cotton pellets within reach during treatment while at the same time protecting them hygienically from spray and dust. A firm base and well-designed lid allow dispensing cotton pellets with only one hand. The rounded rods of the dispensing area prevent the pellets from spueezed on removal as opposed to a dispensing area with a grid structure.

## Order information

	Contents	REF
<b>LUNA</b>	<b>300 g</b>	
	# 1	100 001
	# 2	100 002
	# 3	100 003
	# 4	100 004
	<b>500 pcs.*</b>	
	# 1	105 001
	# 2	105 002
	# 3	105 003
	<b>1000 pcs.*</b>	
	# 1	104 001
	# 2	104 002
	# 3	104 003
<b>LUNA FOR KIDS</b>	<b>200 pcs.</b> # 2	100 032
<b>LUNA GOLD</b>	<b>300 g</b>	
	# 1	109 121
	# 2	109 122
	# 3	109 123
	<b>1000 pcs.*</b>	
	# 1	109 101
<b>COTTON-ROLL</b>	# 2	109 102
	# 3	109 103
	<b>1000 pcs.*</b>	
	# 1	104 011
<b>ROLL-O-MAT</b> plastic, white 134 × 100 × 48 mm	# 2	104 012
	# 3	104 013
	<b>1 pc.</b>	120 004
	<b>chrome nickel steel</b> 134 × 100 × 43,5 mm (can be used with cotton roll sizes 1–3), standard length (36 – 40 mm). Not suitable for LUNA FOR KIDS.	120 005
<b>LUNAMAT</b> plastic, white 52 × 106 × 100 mm (can be used with cotton roll sizes 1–4), standard length (36 – 40 mm) and with LUNA FOR KIDS.	<b>1 pc.</b>	120 000

\*Only for export

© 2016 Coltene/Whaledent AG - www.coltene.com



			Contents	REF
<b>PAROTISROLL</b>	<b>ø</b>	<b>L</b>		
	# 1	9 mm	80 mm	100 pcs.
	# 2	10 mm	80 mm	160 002
	# 3	10 mm	100 mm	160 003
	# 5	8 mm	100 mm	160 005
	anatomic	10 mm	100 mm	50 pcs.
<b>PAROTIMAT</b> 112 × 103 × 82 mm			<b>1 pc.</b>	162 100
<b>COTTON PELLETS</b>	<b>ø</b>			
	3.0 mm	# 000	1.5 g	210 000
	4.0 mm	# 00	4.0 g	210 100
	4.8 mm	# 0	4.0 g	210 110
	6.3 mm	# 1	10.0 g	210 111
	7.6 mm	# 2	10.0 g	210 112
<b>SOLOMAT<sup>1</sup></b> stainless steel 72 × 53 × 50 mm		# 000	1 pc.	6001 9711
		# 00		6001 9712
		# 0		6001 9713
		# 1		6001 9714
		# 2		6001 9715
<b>SOLOMAT-CLICK<sup>2</sup></b> plastic 43 × 50 × 41 mm		# 000	1 pc.	224 000
		# 00		224 100
		# 0		224 110
		# 1		224 111
<b>REFILLS</b>		# 000	3 pcs.	214 000
		# 00		214 100
		# 0		214 110
		# 1		214 111
<b>SOLOMAT-N<sup>2</sup></b> plastic 75 × 43 × 52 mm		# 000	1 pc.	223 000
		# 00		223 100
		# 0		223 110
		# 1		223 111
		# 2		223 112
<b>REFILLS</b>		# 000	6 pcs.	213 000
		# 00		213 100
		# 0		213 110
		# 1		213 111
		# 2		213 112
<b>COMPRECAP</b>		# 1	120 pcs.	530 008
		# 1	600 pcs.	531 008
		# 3	120 pcs.	530 010
		# 3	600 pcs.	531 010
		# 5	60 pcs.	530 014
		# 5	400 pcs.	531 014
		Set (#1/3/5)	120 pcs.	530 000
<b>COMPRECAP ANATOMIC</b>		# 1	120 pcs.	531 201
		# 3	120 pcs.	531 203
		# 5	60 pcs.	531 205
		Set (#1/3/5)	120 pcs.	531 299

<sup>1)</sup> Dispenser is supplied without cotton pellets

<sup>2)</sup> Dispenser is supplied with cotton pellets

### Coltene/Whaledent Ltd.

The President Suite  
Kendal House, Victoria Way  
Burgess Hill, West Sussex  
RH15 9NF / UK  
T +44 1444 235486  
F +44 1444 870640  
info.uk@coltene.com

### Coltene/Whaledent

GmbH + Co. KG  
Raiffeisenstraße 30  
89129 Langenau / Germany  
T +49 7345 805 0  
F +49 7345 805 201  
info.de@coltene.com

### Coltene/Whaledent Inc.

235 Ascot Parkway  
Cuyahoga Falls, Ohio 44223 / USA  
T +1 330 916 8800  
F +1 330 916 7077  
info.us@coltene.com

